

## Instructions for Removing the Water Flosser Water Tank:

### 1. Power Off the Device:

- Ensure the water flosser is turned off.


### 2. Locate the Water Tank:

- Find the water tank, typically located at the back or bottom of the water flosser unit. It's usually transparent and has a fill cap.



### 3. Release the Water Tank:

- Holding the unit with the front part facing away from person, pull the water tank downwards towards the front of the main unit.

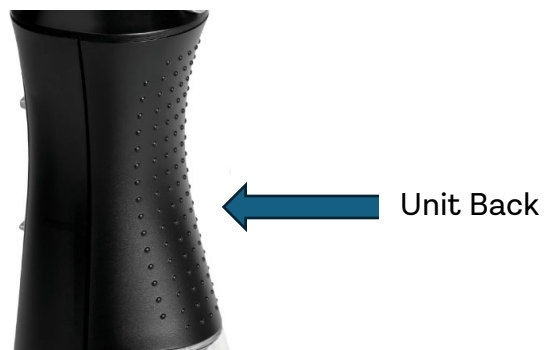

### 4. Remove the Water Tank:

- Once released, carefully lift or slide the water tank off the unit.

**5. Empty the Tank (if needed):**

- If the tank is full or has leftover water, pour it out before cleaning or refilling.

**6. Reinstall the Water Tank:**

- After cleaning or refilling, align the tank back onto the device, ensuring it fits securely. You should hear or feel it click into place.

**Tips:**

- Be gentle when removing the water tank to avoid damaging any connecting parts.
- Some tanks may have rubber gaskets or seals. Ensure these are properly aligned when reattaching the tank.
- Clean the tank regularly to prevent the buildup of bacteria and mineral deposits.