



1 Explanation of symbols

	See chapter 2. Safety
	Important information for the user
	Product can be sterilized in a steam sterilizer (Autoclave)
	Can be used in thermal disinfectors
	CE-mark – Confirms that this product fulfills all requirements for medical products
REF	Reference number
SN	Serial number

2 Safety

2.1 Description of danger levels



	CAUTION
	CAUTION is used if a deficiency in adequate care may result in an endangering of the patient, the user or others.

	WARNING
	WARNING is used if improper use may result in serious injuries.

2.2 Safety notes

	WARNING
	Risk of injury for the patient and the operator Concerning damages like irregular running noise, irregular vibrations, an unspecific rise in temperature and/or other defects. → Stop operating immediately and contact our service.

	WARNING
	Risk of infection In case of an injury of the tissue in the oral area, do not proceed to work with air-operated instruments because of a high risk of an infection.

 WARNING	WARNING
	Burns due to an increased heat of the head of the instrument Burns in the oral area may occur due to an increased heat of the head of the instrument.
 CAUTION	CAUTION
	Accelerated deterioration and dysfunctions due to long term non-usage → Before storage, clean, maintain and dry the instrument according to the norm. Store in the original packing, dry and dustproof.

3 Product description

This manual can be used for the following BENCO products:

"e-type Airmotor Handpiece Swing Latch Head" [4104-345]	
---	--

3.1 Intended purpose


This medical product may only be used in combination with the suitable contra angle body and for the treatment of patients in dental surgeries.

3.2 Use according to regulations


This product may only be used for the intended purposes, executed by professionally trained personnel. The indications and warnings in this manual must be considered. This also includes the observing of the processing (see chapter **5 Processing according to ISO 17664**) and the valid national health and safety regulations. Any other use may result in an endangering of the patient, the operator or others and is strictly prohibited. Before each application, the perfect functionality of the instrument must be tested.

3.3 Technical specifications

"e-type Airmotor Handpiece Swing Latch Head"	[4104-345]	
Transmission/reduction	1:1	
Max. speed	20 000 min ⁻¹	
Suitable bodies	Any 12 tooth (NSK Style), contra angle bodies	
Suitable burs/attachments	Contra angle burs	
Manual Latch	yes	

 NSK® is a registered trademark. The Manufacturer does not have any economical connection to the company mentioned above.

3.4 Transport and storage requirements



 CAUTION	CAUTION
	Dysfunction of the instrument while operating after an overcooled storage → Heavily cooled instruments must be warmed to room temperature prior to operating.


For transportation and storage, please consider the following:

Temperature	-20 °C - +70 °C (-4 °F - +158 °F)
Relative humidity	5% - 95%, non-condensing
Air pressure	700 hPa - 1060 hPa (10 psi - 15 psi)


Protect this medical instrument from moisture!


4 Putting into service and operation

 WARNING	WARNING
	Risk of infection for patients and operator due to non-sterile instruments → Before the first and after each further use the instrument must be processed and sterilized as defined in chapter 5 Processing according to ISO 17664 .
 WARNING	WARNING
	Risk of infections for others at waste disposal → Before waste disposal, this medical product has to be processed according to instructions and must at least be sterilized.

 CAUTION	CAUTION
	Dysfunctions and/or risk of an infection caused by moist and contaminated compressed air and by contaminated cooling water → This product may only be operated with CE-certified dental units or with others certified by actual national guidelines and with a corresponding supply of water and compressed air according to the valid norms of dental care.


4.1 Attaching and removing the contra angle head


 WARNING	WARNING
	Risk of injury for the patient due to a separating of the contra angle head An incorrectly attached head may loosen while operating. → Prior to each treatment check if the head is firmly attached and if the coupling nut is completely tightened.

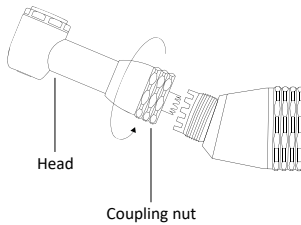



- ☞ Hold the head as shown in the figure and insert it.
- ☞ Tighten up the coupling nut in direction of the arrow.
- ☞ To remove the head unscrew the coupling nut in opposite direction and take off the head.

4.2 Insertion and removal of rotating instruments (LBH11)

 WARNING	WARNING
	Risk of injury for the patient due to a falling out of the rotating instrument → After each insertion of burs check the holding force of the spindle by pulling on the bur.


 WARNING	WARNING
	Risk of injury due to the use of worn out rotating instruments Rotating instruments can fall out while operating and thus injure the patient. → Do not use any rotating instruments with a worn out shaft.

 WARNING	WARNING
	Risk of injury for the user due to rotating instruments Cuts and resulting infections may occur. → Always use gloves or a finger protection.

 WARNING	WARNING
	Danger due to a defective spindle In case of a defective spindle, the bur may fall out and cause injuries. → Check the correct attaching of the spindle by slightly pulling at the bur. → Use gloves or finger protection for insertions and removals of the spindle to prevent injuries and infections.

i Only use rotating instruments made of steel or carbide metal, that fulfill the EN ISO 1797-1 requirements and also fulfill these criteria:


Diameter of the shaft	2.334 to 2.350 mm
Full length of the shaft	max. 22 mm

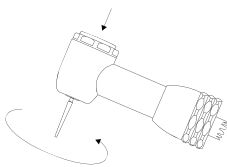

- ☞ Move the swing latch outward to the open positions and insert the rotating instrument into the spindle with turning movements.
- ☞ Afterwards release the move the swing latch to the lock position.
- ☞ Check the proper fitting by slightly pulling at the bur.
- ☞ To remove the head let the bur come to a full stop. Move the swing latch to the open position bur can be released and then be removed.


5 Processing according to ISO 17664

5.1 Processing at the location of use


 WARNING	WARNING
	Risk of infection due to non-sterile instruments. → Always wear gloves.


- ☞ Process the instrument within one hour after each treatment in order to prevent piling.
- ☞ Immediately remove any debris of cement, blood or composite.
- ☞ Ensure a dry transport of the instrument to the place of processing.
- ☞ Do not put the instrument in any kind of solution.

5.2 Cleaning

 WARNING	WARNING
	Dysfunctions due to cleanings in an ultrasonic bath → Only clean the instrument as described below.

5.2.1 Manual pre-cleaning

 CAUTION	CAUTION
	Risk of infection due to multiple uses of dental brushes → To prevent cross infections always use a new dental brush for each instrument.


 CAUTION	CAUTION
	Insufficient cleaning of the surface and the water- and air pipes → In order to ensure a successful cleaning a thorough manual pre-cleaning, even with the small brush, must be done.

Tools:

- Drinkable water 30 °C ± 5 °C (86 °F ± 10 °F)
- Brush or soft tooth brush
- Dental brush
- 50 ml syringe without cannula

- ☞ To pre-clean the instrument place it under running water and carefully brush debris completely off the surface by using one of the tools mentioned above.
- ☞ To pre-clean the spindle, please insert the dental brush in the spindle opening. Carefully remove any residue and debris by circular motions with the brush.
- ☞ Repeat the process several times and each time clean the dental brush with water.
- ☞ To pre-clean the water- and air pipes place the syringe at the lower end of the handpiece and rinse the cavities at least 5 times.

5.3 Disinfection

 CAUTION	CAUTION
	The use of disinfecting baths and/or of chlorinated disinfectants may lead to defects and dysfunctions of the instrument. → Only clean and disinfect the instrument in a thermal disinfectant.

5.3.1 Manual disinfection of the surface

Only use disinfectants whose microbiological effectiveness is ensured by the manufacturer (e.g. VAH/DGHM-registration and CE-mark). Concerning the material compatibility, only disinfectants free of aldehyde and based on alcohol may be used.

- ☞ Apply the disinfectant to a cloth and wipe the instrument clean.
- ☞ Consider the reacting time that is set by the manufacturer of the disinfectant.
- ☞ Follow the manual for the disinfectant.

5.3.2 Mechanical cleaning and disinfection of the interior and the exterior



Benco recommends thermal disinfectors that fulfill the EN ISO 15883-1 requirements and are approved by the manufacturer for the straight- and contra angle, as well as for high speed handpieces. The cleaning should be carried out at a minimum of 55°C (131°F) and for at least 5 minutes, the disinfection at a minimum of 90°C (194°F) and for at least 5 minutes (for an A0-value > 3,000). For the cleaning a mild alkaline cleaning agent with a pH of 9 to 11 is recommended, like e.g. Neodisher® MediClean Forte. For the validation of the process see chapter 6.2 **Validation of the processing.**

- ☞ For the internal cleaning the adaptors specified by the manufacturer must be used.
- ☞ Repeat the process, if visible contaminations still exist after the treatment in a thermal disinfectant.
- ☞ Before use the instrument must be dry and free of residues.
- ☞ To prevent any kind of impairments of the instrument, make sure after every cycle that the instrument is dry on the inside and the outside.
- ☞ Directly afterwards use the Vector Premium Service Oil [4434-533] to lubricate the instrument.

5.4 Process of drying

5.4.1 Manual drying

☞ Blow on the instrument with compressed air from the inside and the outside until no visible drops are left.


5.4.2 Mechanical drying

The drying process normally is part of the cleaning process of your thermal disinfectant. Please follow the instructions for your thermal disinfectant.

☞ In case of moisture remaining on the instrument after the cleaning process, blow with compressed air until completely dry.

5.5 Maintenance with the Vector Premium Service Oil VL-S [4434-533]

Lubricate and maintain the instrument after each mechanical cleaning and before each sterilization process.

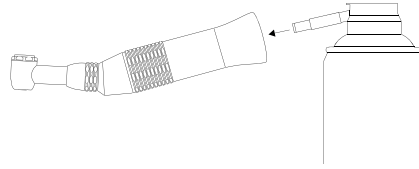

☞ For the maintenance of the cartridge screw the enclosed lubrication tool onto the lubrication can.

☞ Insert the lubrication tool into the lower end of the instrument and spray into the canal for drive air for approximately 2 seconds until the oil evades from the head of the contra angle handpiece.



☞ Plug the instrument onto your treatment unit and let it run for 15 to 20 seconds so that the surplus lubricant can escape from the instrument.


6 Sterilization

i The sterilization bag has to be big enough for the instrument so that the instrument fits without stretching the bag. Concerning quality and application, the bag has to fulfill all the valid norms and has to be suitable for the sterilization process.

☞ Shrink-wrap each instrument individually into a sterilization bag.

6.1 Sterilization with a steam sterilizer (Autoclave) according to EN 13060/ISO 17665-1

 CAUTION	CAUTION
	Damages due to inappropriate maintenance or service may result in dysfunctions and/or deterioration. Reduced lifetime of the product → Lubricate the instrument with the Vector Premium Service Oil [4434-533] before each sterilization.
 CAUTION	CAUTION
	Contact corrosion caused by moisture Damages of the product → After sterilization, immediately remove the instrument from the autoclave.

 This medical product has a temperature resistance of up to max. 135°C (275°F).

Only apply the following sterilization process:

- Autoclave with a triple pre-vacuum (recommended): minimum of 5 minutes at a minimum of 134°C (273°F)

☞ Pay close attention to the manufacturer's instruction.

6.2 Validation of the processing

i The validation of the cleaning and of the disinfection was carried out with the device Miele® G 7835 in the program Vario TD with the cleanser Neodisher® MediClean Forte.

Adjustment:

- Pre-rinsing for 1 minute

- Cleaning at 55°C (131°F) for 5 minutes with a dosage of 0.5% Neodisher® MediClean Forte


- Neutralization

- Washing down

- Disinfection at 90°C (194°F) for 5 minutes

The sterilization was validated with the device Tuttnauer® EHS 3870 in a fractional pre-vacuum process at 134°C (273°F) for 5 minutes. Please note, that the reprocessing procedure performed in your dental office must be validated also.

6.3 Maintenance & service

 CAUTION	CAUTION
	Inappropriate maintenance or service may result in dysfunctions and/or deterioration. → Always perform proper maintenance and service.

7 Storage

- ☞ Store the cleaned and sterilized instruments in a dry, dark and cool place, sheltered from dust and germs.
- ☞ Please pay attention to the date of expiry of the sterilization liquid.

8 Tools & spare parts

Contra angle head	[4104-345]	
Lubrication tool		VL-AS
Premium Service Oil		VL-S [4434-533]
Premium Service Oil for KaVo® QUATTROCARE		VL-AQ [4374-778]



KaVo® is a registered trademark. The manufacturer does not have any economical connection to the company mentioned above.

9 Warranties

Warranties for this product are stated in the Terms of Purchase on Benco.com.

Distributed by

BencoDental™

800.GO.BENCO • BENCO.COM
295 CenterPoint Blvd, Pittston, PA 18640

Made in Germany