QuickPass(? IN-OFFICE DENTAL WATER TEST

Routine testing of dental unit waterlines is the only way to know if you're delivering safe water to patients. Not only is it recommended by the CDC, ADA, & OSAP, but because 31% of treated waterlines still fail to meet the <500 CFU/mL standard, testing is key to protecting your patients and practice[†]. The brand new QuickPass™ In-Office Dental Water Test is more, and easier to read than any other in-office dental water test.





The only in-office water test designed specifically for dentistry



Most reliable media with germicidal neutralization included



Fastest results with 48 – 72 hour incubation period



Advanced media makes bacteria grow red for easier reading



The most cost-effective waterline testing solution available



888.843.3343 | proedgedental.com







Study cited in "Treating & Monitoring Dental Water", Compendium of Continuing Education in Dentistry, by Nancy Dewhirst, RDH, BS & John A. Molinari, PhD. Data collected from 22,196 consecutive waterline tests (R2A, neterotrophic plate count tests) conducted in 2017 by ProEdge Dental Water Labs. All product types provided by dental offices. Averages may not reflect your test results. All Pass/Fail measured against CDC's <500 CFU/ml standard



1.800.GO.BENCO or visit bencodental.com



How often does my practice need to test?


The Organization for Safety, Asepsis, and Prevention (OSAP) makes the following recommendations for testing frequency of dental unit water quality:

- If available, follow your dental unit and treatment product manufacturer recommendations.
- > Otherwise, test waterlines every month until passing results are consistent for two consecutive months, then test at least quarterly.
- > If your waterlines fail a test, take corrective action and retest to verify remediation was effective.
- > Test following any changes to your treatment protocol, equipment, source water, or plumbing.

Quarterly Waterline Testing Demonstrates CDC Compliance to Protect Your Practice:

Documentation of guarterly dental water testing proves an effective treatment protocol and your office's compliance with the CDC guideline for dental unit water quality of <500 CFU/mL. With passing results, you have protection in the event of an inspection, infection, or complaint.

Sample QuickPass Testing Log. Download at ProEdgeDental.com/QP-Log

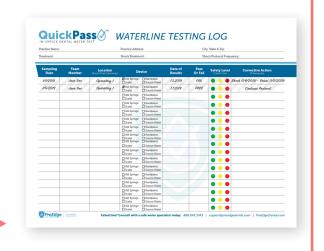

Questions?

888.843.3343 | ProEdgeDental.com | support@proedgedental.com



YOUR ASSURANCE