



Obturation Success for Everyone



Introducing BioRoot® Flow Bioceramic Root Canal Sealer



Easy

Ready-to-use syringe provides easy extrusion directly into canals and you maintain your preferred obturation technique.

Successful

Highly-pure calcium silicate mineral formulation (resin-free) means no shrinkage⁽¹⁾, limited bacterial growth, and increased biocompatibility and bioactivity which reduce the risk of adverse reactions.

For Everyone

Designed for cold and warm techniques with excellent flow into accessory canals⁽²⁾, tight adhesion to dentin, and gutta-percha, and reproducible results in all root canal treatments from simple to complex.⁽³⁾

Septodont

205 Granite Run Dr., Suite 150
Lancaster, PA 17601

800.872.8305

www.septodontusa.com

(1) Internal data. Pr Colon report

(2) Drukteinis & Camilleri : Bioceramic materials in clinical endodontics. Springer 2020

(3) BioRoot™ flow SEM interfaces. Dr Camilleri


BencoDental™



Technical Insights

Innovative tip allows direct and precise placement

- ▶ Flexible tips allow bending and access to all root anatomies
- ▶ Precise 21-gauge diameter for optimal delivery limiting the risk of going too deep

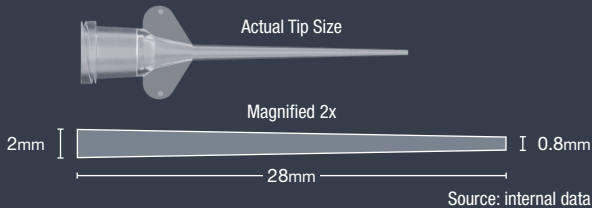


Highly Radiopaque

- ▶ 6.3mm Al radiopacity
- ▶ Easily visible on X-rays

Easily Retreatable

- ▶ Easy and fast to remove if necessary
- ▶ Use hand files or rotary NiTi Files
- ▶ BioRoot Flow is removed to working length in less than 10 minutes.

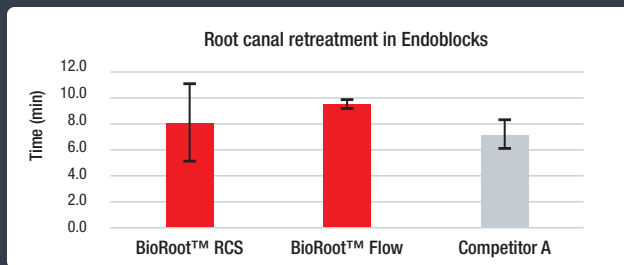

Source: internal data

Septodont


205 Granite Run Dr., Suite 150
Lancaster, PA 17601

800.872.8305

www.septodontusa.com

Long Lasting pH

- ▶ High pH is maintained over time creating an alkaline environment in which bacteria will not grow

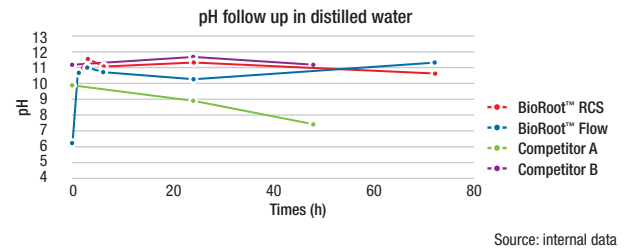

Source: internal data

Technical Features

- ▶ Working time | > 60 min
- ▶ Setting time | 5 to 6 hr

Source: internal data; Dr Camilleri

Product Information

- ▶ (1) 2g syringe
- ▶ (20) intra-oral tips
- ▶ (1) finger grip
- ▶ Item #5962-005

